# SN J6802B-02A

# S E R V I C E N O T E

Supersedes: J6802B-02

## J6802B DNA MX

Serial Numbers: All units prior to SN SG44320212 for Option 001 SG44330157 for Option 002

Line filter improperly wired in relation to fuses.

To Be Performed By: Agilent-Qualified Personnel or Customer

**Parts Required:** 

P/N Description Qty.

2110-0003 or equiv. Fuse 3A 250 V

### ADMINISTRATIVE INFORMATION

SERVICE NOTE CLASSIFICATION:			
MODIFICATION RECOMMENDED			
ACTION CATEGORY:	[[]] IMMEDIATELY [[]] ON SPECIFIED FAILURE	STANDARDS: LABOR: 0.25 Hours	
	X AGREEABLE TIME		
LOCATION CATEGORY:	X CUSTOMER INSTALLABLE [[]] ON-SITE X SERVICE CENTER	SERVICE [[]] RETURN INVENTORY: [[]] SCRAP X SEE TEXT	USED [[]] RETURN PARTS: [[]] SCRAP X SEE TEXT
AVAILABILITY:	PRODUCT'S SUPPORT LIFE	AGILENT RESPONSIBLE UNTIL: Dec 2007	
AUTHOR: DM	PRODUCT LINE: 2J		
ADDITIONAL INFORMATION:			

1

© AGILENT TECHNOLOGIES, INC. 2006 PRINTED IN U.S.A.



Page 2 of 6 J6802B-02A

#### **Situation:**

Do to an error in the assembly documentation the line filter was incorrectly wired on most of the first 150 units shipped. The fuse provides no protection for surges in AC line voltages to the power supply. This causes the power supply internal fuse to open and then the power supply ceases to function.

#### **Solution/Action:**

Disconnect power cord.

Using a small thin blade screwdriver, open the fuse-block door to access the fuse block. (See photo #1)

Remove fuse-block from line filter. (see photo #2)

Remove shunt from fuse-block and replace with fuse. (See photos #3 & #4)

Replace fuse-block into line filter and close door.

Reconnect power cord

Page 3 of 6 J6802B-02A

Photo #1 Note: If your unit is EC 11:00 AND has the 115V visible in the window, stop here as this unit already has two fuses.



Page 4 of 6 J6802B-02A


Photo #2



Slide fuse-block (red plastic) back towards you. You may need to use the screw driver again.

Page 5 of 6 J6802B-02A

Photo #3

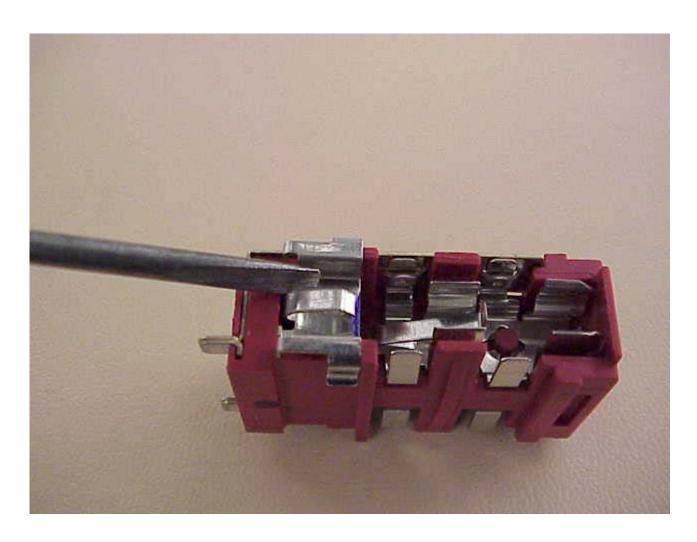

Page 6 of 6 J6802B-02A

Photo #4

